


User Manual







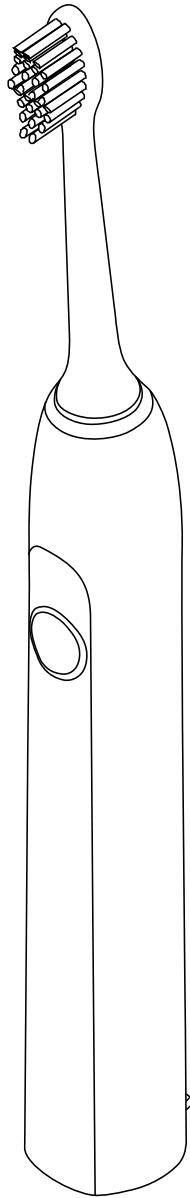





**CE, FCC, FDA, RoHS, EMC,
IPX7**

TO OUR VALUED CUSTOMERS:

First of all, thank you for choosing our sonic electric toothbrush, which affords us an opportunity to provide you with the best quality products and services. This product is a precision intelligent electric appliance, which can effectively remove most of dental plaque and dental calculus. Practice has proved that continuous use can effectively reduce periodontitis, dental plaque and other oral problems so as to make your smile brighter and a good life simpler!

In order to obtain the best experience in using, please read this introduction carefully.

Product features:

	Ergonomic Design		Dynamic Fluid Cleaning Tech.
	Multiple Tooth Protection		Still Strong Under Pressure
	Pressure Sensor		100-Day Battery Life
	Remove 90% Of Stains		Washable
	3 Levels Charge Display		Easy To Carry

WARNING

- Please make sure the charger is put in dry area, away from the water. Do not use this product while taking a bath or shower.
- Please do not touch the charger with water entering or the one fallen into the water. In such case, disconnect the charger immediately.
- Do not use the charger with a damaged wire or other problems.
- Please do not use this product or its charger out of doors or in high temperature areas.
- Please make sure the voltage of the charger is the same as its outlet. If it is not easy to insert the outlet, do not make it by force.
- If the charger will not be used, please disconnect the charger so as to avoid firing or other dangers.
- Do not replace the battery and power line on your own. If the power line is damaged, please scrap the charger.
- This product does not contain the user-serviceable parts. If the product is damaged or needs to be repaired, please stop using them including the brush head, handle or charger and contact the manufacturer for warranty.
- This product is intended for cleaning teeth, gums and tongues. Do not use it for other purposes. Please use this product under the advice of the dentist.
- This product is designed for people aged eight or above, people (including children) who are physically, sensuously or mentally handicapped or lacking of relevant experience and knowledge need to be supervised and directed to ensure their safe use. Do not let children clean and maintain the product.
- Please keep the product and its cable out of the reach of the children aged below eight years old, take care of the children to ensure that they don't play with the product.
- This product is a personal care device and should not be shared by multiple people in dental clinics or public places.
- If you have undergone surgery on your mouth in two months, please consult a dentist before use.
- If you feel uncomfortable or painful after using this product, with excessive bleeding or a slight one in the gums continuing for more than 1-2 weeks, please consult a dentist.
- If you have a pacemaker or other implanted devices, please consult a dentist before use.


ELECTROMAGNETIC FIELDS(EMF)

This product complies with the safety standards for the electromagnetic devices.

PRODUCT SPECIFICATION

Description: D7 series sonic electric toothbrush
 Category: sonic electric toothbrush
 Model no. D7
 Rated voltage: DC 3.7V
 Rated power: 1.8w
Charger
 Rated voltage: DC 5V
 Rated power: 2.5w
 Execution standards: JSSANQT2017
 Origin: Made In China

PART LIST

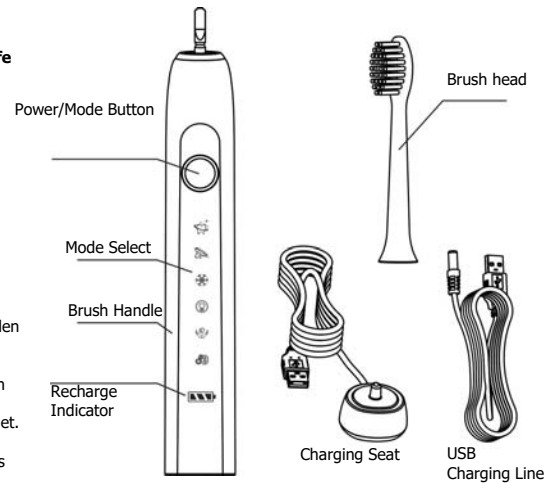

PRODUCT SPECIFICATION

Brush handle	Brush head	Charging seat/Charging line
1	2	1

*The charger way and the accessory vary with the product model.

CHARGING

Charging Seat (for specific model) :

Put the electric toothbrush into the charging seat, plug one end of USB line into the USB charger, the electric toothbrush will be fully charged after 24 hours when the toothbrush indicator turns from red to three-level green (Single-level green for standard products).




USB Charging line with large current (for specific model) :

Insert the small end of the USB charging line into the bottom of the toothbrush, insert the other end of the line into the USB charger, the toothbrush can be fully charged in about 6 hours when the indicator of the toothbrush turns from red to green. 1 minute charging can meet the need for twice use with quick charging.

Recharge Indicator:

The battery charging indicator will show the remaining amount of electricity in the battery

Deluxe three-level charging indicator (for specific model)

	Three-level constant green light indicates electricity amount in the battery is 70-100%;
	Two-level constant green light indicates electricity amount in the battery is 40-70%;
	One-level constant green light indicates electricity amount in the battery is 10-40%

One-level constant red light indicates electricity amount in the battery is below 10% and it needs to be charged as soon as possible.
 One-level flashing red light indicates the battery needs to be charged immediately.


Standard charging indicator (for specific model)

Constant green light indicates electricity amount in the battery is 50-100%;

Flashing green light indicates electricity amount in the battery is below 50%;

Constant red light indicates electricity amount in the battery is below 10% and it needs to be charged as soon as possible.

Flashing red light indicates the battery needs to be charged immediately.

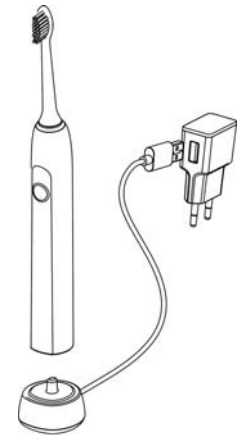

Note:

- If the product is used for the first time or has not been charged for a long time, resulting the toothbrush cannot start and the indicator may not display, which are normal phenomenon, please be patient, the battery needs to be charged for six to ten minutes before the indicator resume display.
- If the toothbrush stops vibrating or doesn't vibrate and the red light of the power indicator flashes, meaning that the electric quantity is insufficient, please charge it timely.
- Please use the charger corresponding to the product model.
- To ensure the electricity amount in the battery is always enough, you may charge the brush before the red light appears.
- It is normal phenomenon that the brush handle is a little bit hot to you; if the handle of the toothbrush feels very hot, please cut off its power immediately and contact the manufacturer.

FUNCTIONS & MODES

In the state of shutdown, the product can start up by the button being pressed shortly.

In 3 seconds after the product starts up, with the indicator of the selected mode highlighting and flashing, the modes of the product can be switched by the button being pressed shortly. After 3 seconds, with the display of the indicator changing from flashing to constant lighting, the product can shut down by the button being pressed shortly.

After the tooth brushing mode selected, the toothbrush will shut down automatically after 2 minutes or 1 minute.

The toothbrush will shut down by the button being pressed long and the tooth brushing mode running last time is set as the first mode. When charging, the toothbrush will stop vibrating and the brushing mode running last time will be set as the first mode.



The standard cleaning mode, cleaning the mouth with the daily common force.



Soft mode, cleaning the mouth with small vibration.



Whitening mode, cleaning the mouth and whitening the teeth with big vibration.



Fast mode, completing the cleaning of the mouth in one minute with big vibration.



Nursing mode, massaging gums.


The toothpaste's buffeting function: Pressing the button for 3 seconds in the running, all mode indicators will be turned off and flash once, this will turn the function on/off. And the toothbrush will output with a weak vibration for the first 3 seconds to prevent the toothpaste from being shaken off before reaching the mouth.

Note:

It is recommended to start brushing from the soft mode in the first two weeks of using this product.

GUIDELINES


Step one: Keep the brush head and the shaft of the brush head in the same direction and fix the brush head vertically on the brush handle with a 1mm interval between the brush head and the brush handle.

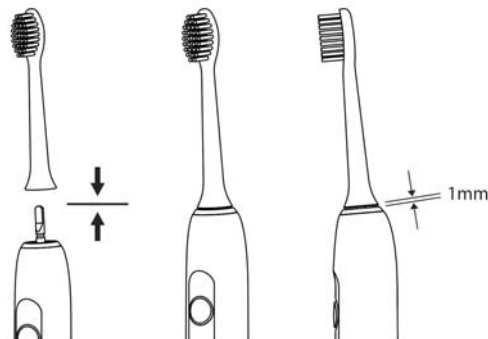


Step two: Wash the head of the toothbrush with water and squeeze a proper amount of toothpaste on the brush head (Put the brush head into your mouth before starting up the toothbrush to prevent the toothpaste splashing).



Step three: Gently place the toothbrush head at the junction of the teeth and gums and make a vertical angle of 45 degrees between toothbrush head and the tooth.

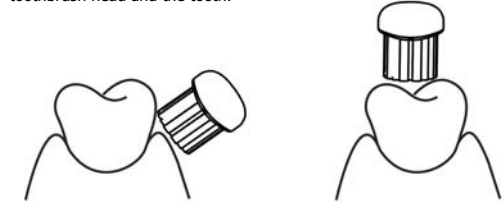

Step four: Press the power button to start up the toothbrush, choose the proper brushing mode, gently hold the toothbrush handle and slowly move the toothbrush along the gum line and the dental occlusal surface.

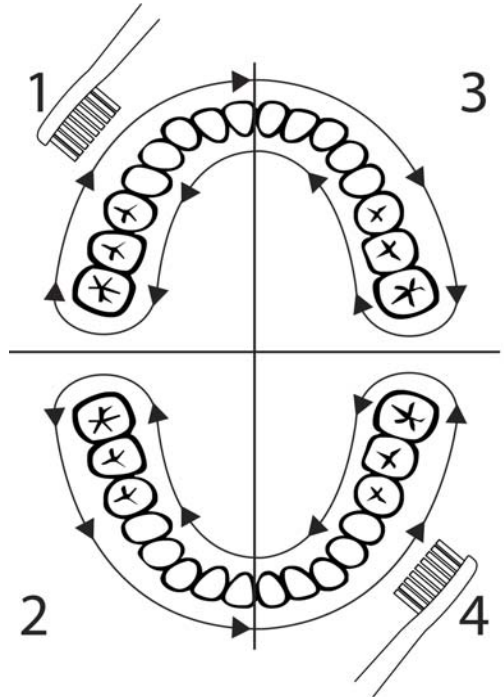

Note:

When in the first time, some people feel tooth irritation or numbness, It is normal, please choose to brush your teeth from soft mode first, and you will adapt to this stimulation within a week. Dentists recommend that you brush your teeth twice a day and not less than 2 minutes each time.

CLEANING & MAINTENANCE

Please clean the brush and brush head after each use and put it in the dry place.



Clean the handle after each brush, and dry it with a soft cloth. Remove the head of the toothbrush from the handle at least once a week, rinse the metal shaft with warm water to ensure that all the toothpaste or residue is cleaned out and dry it with a soft cloth. Do not clean the handle with the water or acid / alkaline cleanser above 40°C. Clean out the stains on the charging seat and adapter with a soft cloth and dry it. Before cleaning the charging seat, make sure that the power plug has been pulled out from the socket, otherwise it might the electric shock.

Note:

- Keep sharp objects away from the rubber seal ring on the metal shaft, otherwise the handle will be damaged.
- Do not clean the toothbrush head, handle, charger and other accessory in the dishwasher, washing machine, etc.
- Do not soak the handle in the water
- Do not push or press the shaft of brush head, which might damage it.

REPLACE THE BRUSH HEAD

The dentist recommends to replace the brush head every three months.

STORAGE

Please disconnect the charger, clean the handle and charger and store them in a cool and dry area.

FAQ

Q: Why doesn't the toothbrush vibrate?

A:

- The electric toothbrush may run out of power, please charge it.
- The socket may disconnect.

Q: Why do I feel itchy when using this sonic toothbrush?

A: The feeling of slight itching / tingling may happen in your initial use, because your mouth has not yet been adapted to this stimulus, usually after using the brush two to three times or for a week, this feeling of discomfort will be relieved.

Q: Why does the toothbrush head vibrate on the teeth?

A: To avoid this problem, please adjust the angle between the brush head and the tooth surface or try to put the brush head on the same plane as the gum line.

Q: Why does the toothbrush sometimes seem to lack of power?

A:


- The toothbrush is running out of power, please charge it;
- The toothbrush head is aging, please replace it.

Q: Does this sonic electric toothbrush require a special charger for traveling?

A: No special charger is needed. The global standard USB interface is used in this electric toothbrush charger, so it can be powered by any power adapter or other digital devices with a USB port.

RECYCLING

A built-in rechargeable battery contains substances that may pollute the environment. Please remove the battery before discarding the product. Please refer to the following way to remove the battery and send the disassembled battery and handle to the recycling centre designated by the local government for environmental protection treatment. If you have problems when you take out the battery, please contact the manufacturer.

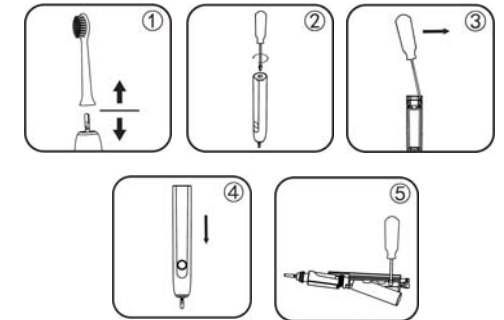

1. Unload the brush head from the handle and turn on the product continuously until it stops vibrating due to battery power is used up.

2. Insert a crosshead screwdriver into the groove located at the bottom of the handle, turn around the screwdriver until the bottom cover open and remove the charging wire rack.

3. Turn the shaft of the brush handle down with its top pressed onto a hard surface, prize off the buckle of the inner component from the edge of the shell with the screwdriver and press down the shell lightly to make the buckle jump out of the groove.

4. Press the shell of the toothbrush downward to make the metal shaft completely pressed into the shell, then take out the inner component.

5. Insert the screwdriver into the plastic bracket of the inner component and pry out the battery.

When following the above steps, be sure to pay attention to the safety operation, protecting your eyes, hands and the work table.

Note:

This operation will irreversibly damage the product.